
# aquadex SmartFlow


## PRODUCT INFORMATION

### AQUADEX CONSOLE

### EASY TO OPERATE, SAFE TO USE

- Adjustable fluid removal rates for customized therapy
- Perform therapy through peripheral or central venous access
- Highly automated with only one setting required to begin
- Hematocrit sensor provides real-time measurement of % blood volume change
- User defined hematocrit limit
- SvO<sub>2</sub> monitoring provides insights into tissue oxygen delivery

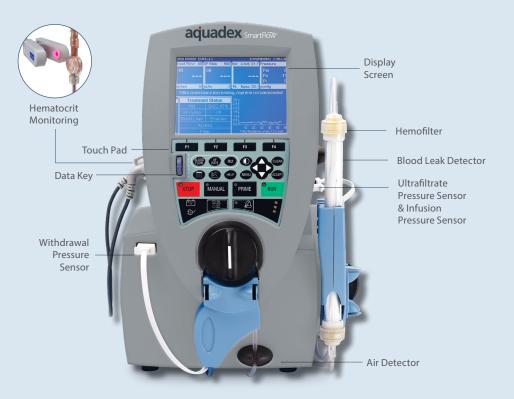

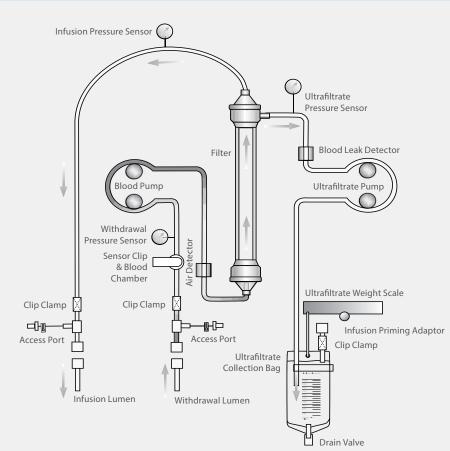

### AQUADEX SYSTEM FLUID PATH

### FULLY INTEGRATED BLOOD CIRCUIT SET

- Patented filter designed to reduce clotting
- Air and blood leak detectors to ensure correct operation
- Needleless access ports for aspiration or infusion of fluids
- Filter Resistance Alert prompts action to help extend circuit life





## PRODUCT INFORMATION

#### **DESIGNED FOR USE WITH DUAL LUMEN EXTENDED THE AQUADEX® SYSTEM** LENGTH CATHETER (DELC) Withdrawal Two independent 16 G lumens Lumen Distal 6 F reverse taper for smooth insertion Opening Infusion Lumen Coil reinforcement for enhanced Distal Opening kink resistance Withdrawal Lumen Clamp (Blue) Withdrawal Lumen Luer Connector Hub Shaft End Caps (Non-Vented) Extension Tubes 15 cm Insertable Length (0.5 cm Depth Markings) Infusion Lumen Luer Connector Approximately 1 cm Infusion Lumen (Not Inserted to Clamp (White) Help Avoid Kinking)

### **ORDERING DETAILS**

MODEL		PRODUCT DESCRIPTION	QUANTITY
English (US) English (UK) French German Portuguese Italian Spanish	120100 120200 120300 120400 120500 120600 120700	Aquadex SmartFlow Console	1
420000		Aquadex SmartFlow Cart/Stand	1
114156 (A2000)		UF 500 Circuit with Hematocrit	1 box of 5

MODEL	PRODUCT DESCRIPTION	SIZE	QUANTITY
114153	Dual Lumen Catheter with Coil Reinforcement	6F (.070" / 1.78 mm l.D.)	1 Box of 5
115463	Catheter with Insertion Kit	6F (.070" / 1.78 mm l.D.)	1 Box of 5 Each

### LEARN HOW AQUADEX WORKS



#### RX ONLY

**INDICATION:** The Aquadex SmartFlow System is indicated for: Continuous ultrafiltration therapy for temporary (up to 8 hours) or extended (longer than 8 hours in patients who require hospitalization) use in adult and pediatric patients weighing 20 kilograms or more whose fluid overload is unresponsive to medical management, including diuretics. All treatments must be administered by a healthcare provider, within an outpatient or inpatient clinical setting, under physician prescription, both of whom having received training in extracorporeal therapies.

Aquadex SmartFlow® is a trademark of Nuwellis, Inc.



Customer Support: 1.855.786.2778 Support@nuwellis.com www.nuwellis.com